

Medical surgical mask differential pressure tester



The picture is for reference only

1.Application

The medical surgical mask pressure difference tester is used to detect the pressure difference of surgical masks. The equipment is equipped with software that can display the pressure curve in real time, output the test report, and save the test data. Meet the test standards YY0469-2011, YY / T0969-2013 ventilation resistance.

2.Features

- 2.1 The equipment is equipped with a pump to simulate artificial respiration.
- 2.2 The detector is equipped with a vacuum sensor and a positive pressure sensor, which are used to collect the test data pressure and reflect it into the computer software.
- 2.3 The pressure has stepless adjustable function, and the technical parameters such as holding pressure and holding time can be set.
- 2.4 Technical parameters of the detector
- 2.5 Flow meter range: 8L / min.
- 2.6 Control method: computer.
- 2.7 Micro-pressure measuring range 0-500pa, accuracy 1pa.
- 2.8 The timing accuracy is 0.1s.





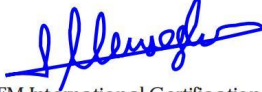
MACHINERY DIRECTIVE ATTESTATION OF CONFORMITY

Technical file of the company mentioned below has been inspected and audit has been completed successfully.

2006/42/EC Machinery Directive Annex VIII has been taken as references for these processes.

Company Name : **JINAN HENSGRAND INSTRUMENT CO., LTD.**
Company Address : NO.4915, WEST JINGSHI ROAD, JINAN, CHINA
Related Directives and Annex : **2006/42/EC Machinery Directive/Annex VIII**
Related Standards : **EN ISO 12100:2010, EN 60204-1:2018**
Product Name : **TESTING MACHINE**
Report No and Date : TCF-JHIC230128-MD
Product Brand/Model/Type : WDS-1, WDS-2, WDS-3, WDS-5, WDS-10, WDS-20, WDS-30, WDS-50, WDS-100, WDS-200, WDS-300, WDS-500, WDW-1, WDW-2, WDW-3, WDW-5, WDW-10, WDW-20, WDW-30, WDW-50, WDW-100, WDW-200, WDW-300, WDW-500, WDW-600, WDW-1000, WE-300, WE-600, WE-1000, WE-2000, WEW-300, WEW-600, WEW-1000, WEW-2000, WAW-300, WAW-600, WAW-1000, WAW-2000, WAW-3000, WAW-5000, WAW-6000, LAW-1000, JB-300B, JBS-300, JBW-300, JB-500B, JBS-500, JBW-500, JBW-750, CST-50, VU-2Y, VU-1S, DWC-40, DWC-60, DWC-80, DWC-100, DWC-196, YES-2000, YES-300, YAW-1000, YAW-300, YAW-2000, YAW-3000, YAW-5000, NDS-200, NDS-500, NDS-1000, NDS-2000, NDW-200, NDW-500, NDW-1000, NDW-2000, NDW-3000, NDW-5000, NDW-10000, HEZ-3, HEZ-6, HEZ-10, JWJ-10, TLS, TLW, HST-P50, HST-P100, HST-P200, HST-P300, RDJ, NDT, DWTT, GBW-60, GW-50, HST, XNR-400

Certificate Number : **M.2023.206.C81366**
Initial Assessment Date : 31.01.2023
Registration Date : 01.02.2023
Reissue Date/No : -
Expiry Date : **31.01.2028**


UDEM International Certification
Auditing Training Centre Industry
and Trade Inc. Co.

The validity of the certificate can be checked through www.udem.com.tr. Upon completion of EC declaration of conformity, it is used solely at the manufacturer's responsibility. This certificate remains the property of UDEM International Certification Auditing Training Centre Industry and Trade Inc. Co. to whom it must be returned upon request. The above named firm must keep a copy of this certificate for 15 years from the registration of certificate. This certificate only covers the product(s) stated above and UDEM must be noticed in case of any changes on the product(s)

Address: Mutlukent Mahallesi 2073 Sokak (Eski 93 Sokak) No:10 Çankaya - Ankara - TURKEY
Phone: +90 0312 443 03 90 **Fax:** +90 0312 443 03 76
E-mail: info@udem.com.tr www.udem.com.tr



HST GROUP Co., LTD.

Jinan Hensgrand Instrument Co., Ltd.

NO. 4915, West Jingshi Road, Jinan 250012, PR CHINA.

Phone/WhatsApp: +86-159 1008 1986

www.hssdgroup.com - e-mail: admin@hssdtest.com

No. 4915, Jingshi Road, Jinan, China

Tel: +86-15910081986